



Interchangeable Docks

The **NEW EZ-Dock Station** integrates four stations into one streamlined unit. It includes a built-in flow circuits and preset timers specifically designed for each cage type, ensuring full AVMA compliance. Each station is independent of each other. This allows for easy docking and removal without having to wait.

Features:

- Compliant with AVMA Guidelines
- Cycle starts automatically when the cage is docked
- Select between adult and neonate settings
- LED indicator lights show cycle status (running/completed)
- Automatically shuts off gas if the cage is removed
- Uses house or bottled CO₂ supply (15-100 psi)
- Docks are interchangeable, offering full customization
- Each station is independently controlled



The dock station can stack up to three rows high, allowing for simultaneous docking of up to 12 cages. The docks are also interchangeable, accommodating future cage model changes. It is currently compatible with specific cage types, including Allentown Nexgen, Nexgen Max, Tecniplast Greenline and Tecniplast Emerald.




Each station is independently controlled




Fits our E-25000
Mobile Workstation

Input: 120-240VAC
1Amp Max @120VAC (per level)
Max Pressure:100 PSI

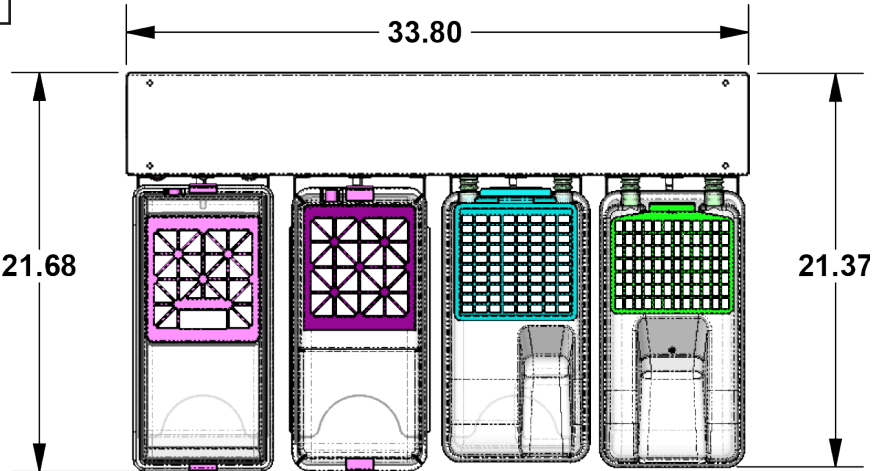

TOP VIEW

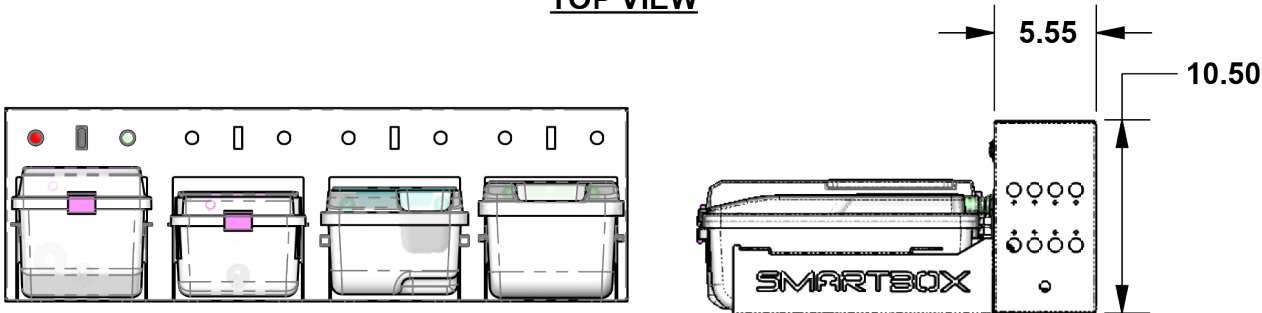

FRONT VIEW